

Silhouette[®]
INTERNATIONAL

Fusion

Glazing | Fusion

Optimal lens base curve is 5.0 Other curves are possible and will require frame adjustments


Glaze with Silhouette mounting pliers with blue tip and BLS 58/98/97



BLS 58 STANDARD	BLS 98 „long“	BLS 97 „short“
Standard Glazing 	for very thick lenses 	for very thin lenses
lens size 2,0 - 6,5 mm	6 - 9,5 mm	1,3 - 3 mm
∅ 1,4 mm	∅ 1,4 mm	∅ 1,4 mm

Inclination | Fusion

1. Open temple
2. Support lens fixation firmly with mounting pliers
3. Hold temple at hinge and adjust inclination at the end piece to the desired position

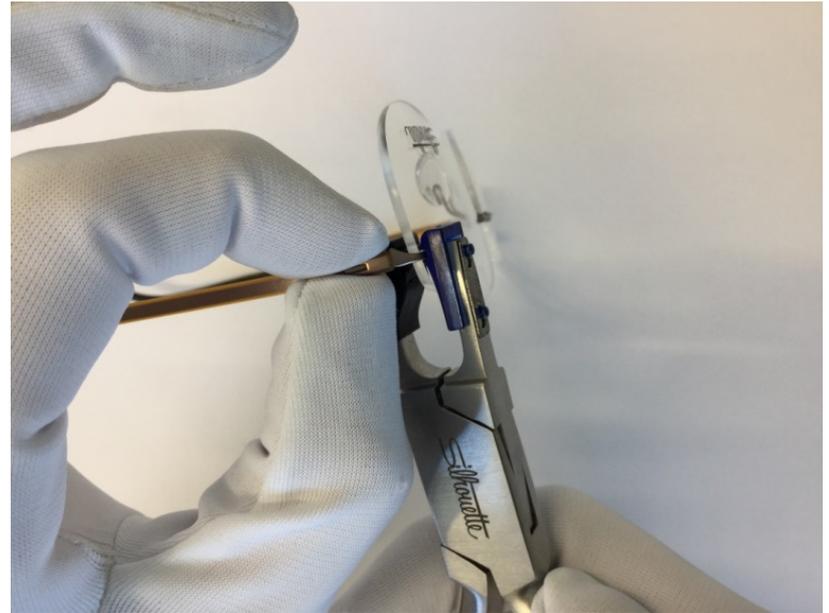

Temple Opening | Fusion



Reduce



Enlarge



Temple disassembly | Fusion

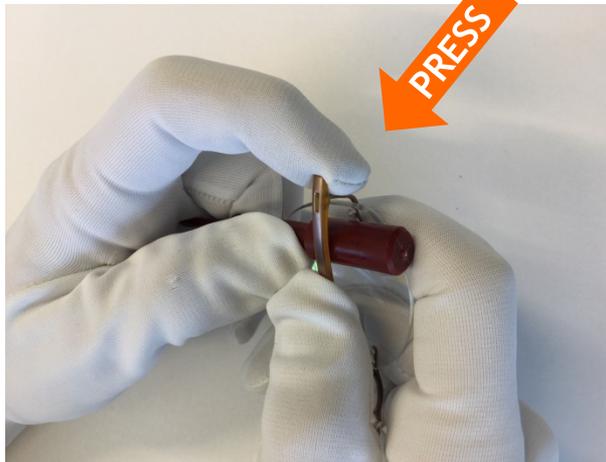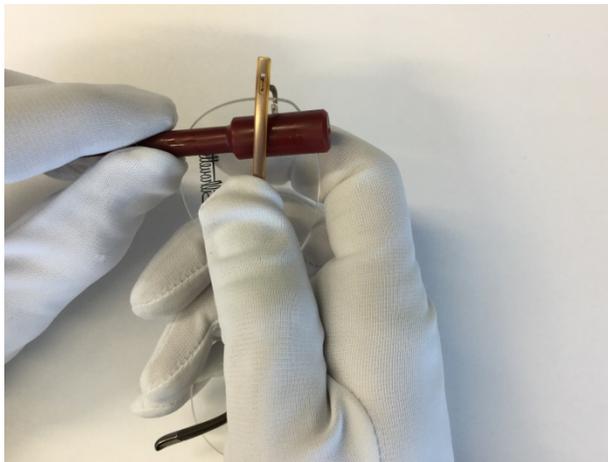

1. Place red multi-tool at the hinge as pictured. Select suitable diameter of multi-tool, according to distance between temple and lens
2. For best protection, please place a cloth between red multi-tool and lens
3. Close temple and apply pressure behind round part.
4. Temple will loosen with a click

Temple Assembly | Fusion

1. Place end piece exactly at holes of SPX-temple
2. Apply counter pressure to endpiece
3. Press SPX-temple into hinge
4. Press the temple inward as shown in picture



Adjust Temple End | Fusion

Metal core inside temple cover allows cold adjustment of temple bend

Bend temple cover in large curves

