



# TECHNICAL INFORMATION (T007-T009; n03)

ROSA – T007  
GEORG – T008  
LINDA – T009  
Spare Part Group n03

# GLAZING

## RECOMMENDATION FOR OPTICAL LENSES:

- Frames have 5.0 base curve
- Flatter lenses should have 5.0 curve bevel applied to lens edge to maintain proper frame curve and ease adjustment

## GLAZING AID:

- Patterns are available through customer service
- OMA shape files are available at [www.silhouettelab.com](http://www.silhouettelab.com)

**Attention: COLD GLAZING ONLY**



# INCLINATION




**Incline the temple at the hinge!**


Special radius design to prevent gap between temple and front when adjusting for pantoscopic tilt.



# TEMPLE ADJUSTMENT

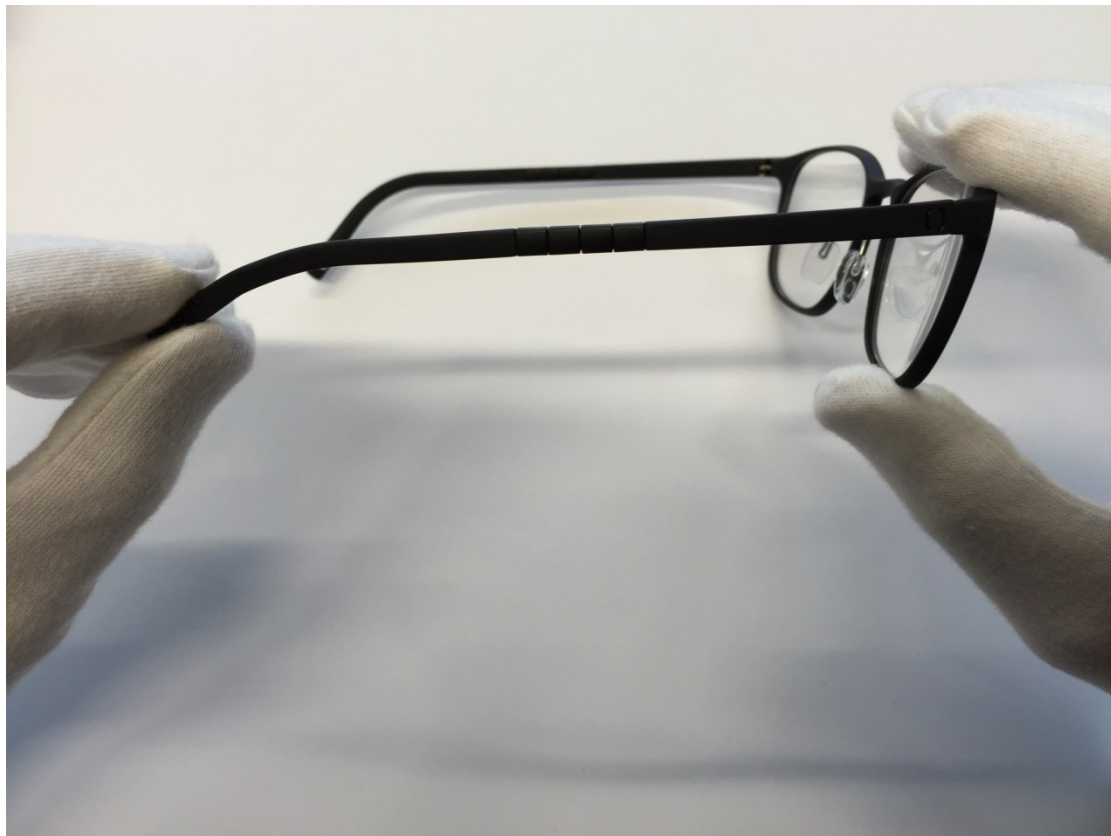

Straighten the temple using small steps

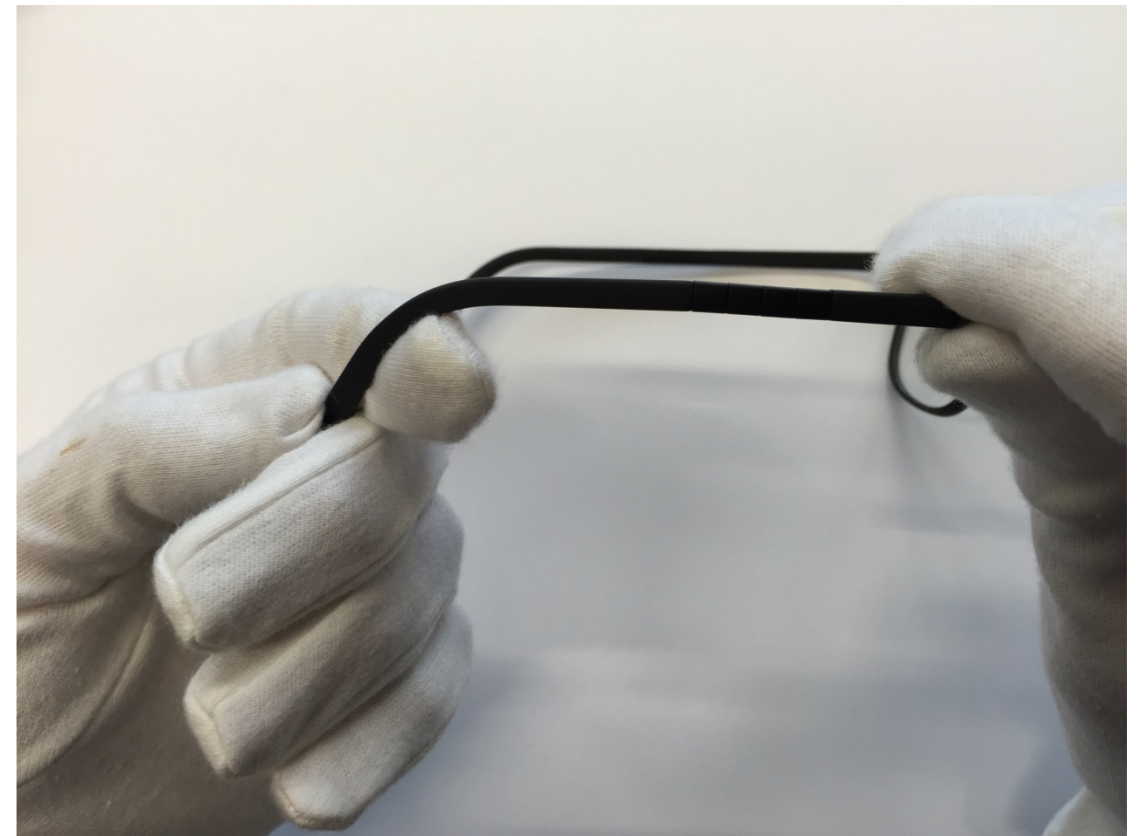

Adjust the temple end by bending in a gentle curve

You may heat the material up to 158 degrees Fahrenheit to ease adjustment